


Device and controls



- 1 Red light
- 2 RF and EMS head
- 3 EMS level indicator
- 4 EMS level button
- 5 Mode indicator
- 6 Power / Mode button
- 7 Charging indicator
- 8 Connector for charging stand
- 9 USB-C charging port
- 10 Charging port on the stand
- 11 Charging stand

Explanation of symbols

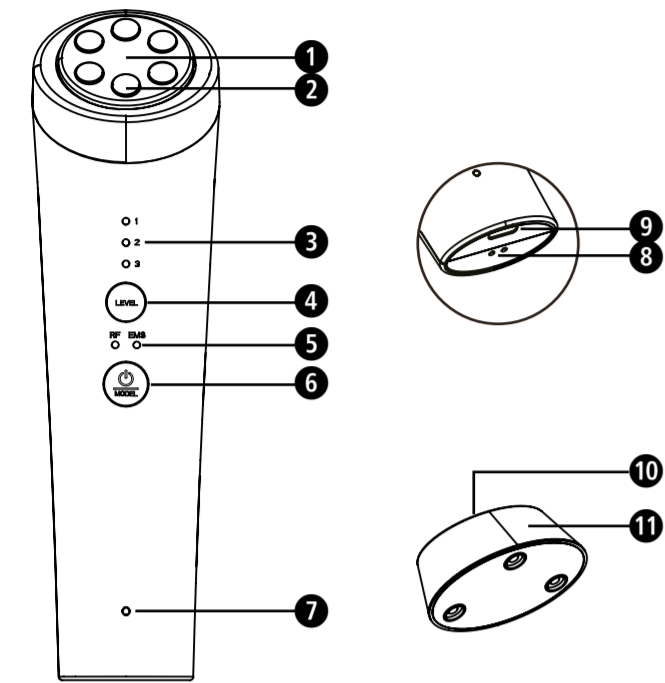

IMPORTANT!
Non-observance of these instructions can result in serious injury or damage to the device.



WARNING
These warning notes must be observed to prevent any injury to the user.



CAUTION
These notes must be observed to prevent any damage to the device.



NOTE
These notes provide useful additional information on the installation or operation.



LOT number



Manufacturer



Recycling symbols/codes:
These are used to provide information about the material and its proper use and recycling.

GB Safety instructions



Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.



- Users of the following medical devices are not allowed to use this device:
 - Users of electronic medical devices that are worn inside the body such as cardiac pacemaker.
 - Users of electronic medical devices for sustaining of life such as artificial cardiopulmonary machine, etc.
 - Users of electronic medical devices attached to the body such as electrocardiogram (ECG) device, etc.
- Patients with the following symptoms or those corresponding to the following should resume using the product after consultation with medical specialist.
 - Those with skin rash or experiencing excessive stinging experience.
 - Those with inflammation as the result of suppuration of the site of acne.
 - Those with particularly sensitive skin with allergic physical constitution.
 - Pregnant and menstruating women, and those who are being treated as an outpatient or has health problems.
 - Those taking hormonal medicine. Those with sensitive skin may experience mild skin troubles upon the use of the device. If such symptom arises, stop using the device for 2-3 days before using the product again.
- Do not use directly on eyes.
- If your skin blisters or burns, stop using the device immediately.
- Make sure to take off all jewelry from the body before using the device.
- Do not store or use in places with high temperature and humidity such as sauna and bathroom, etc.
- Stop using the device if you feel any abnormality in the device during the usage and turn off the power. You may make inquiries to the store where you purchased the product.
- Take precaution as there may be situations of erroneous operation of this device when it comes in contact with other electrical products.
- Stop using the device if there is thunder and lightning.
- Do not use while charging.
- This appliance can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Do not overexpose to the sun after use so as not to cause secondary damage to the skin.
- The following persons shall not use this device: Sunburn patients, healthcare workers with diseases involving photosensitivity and people receiving photosensitive medications.
- Unexpected side effects, such as itching, within 48 hours of first use of the appliance should seek medical advice before further exposure.

About the power supply

- Only charge the unit's battery using the USB charging cable supplied. Never leave the unit unattended when it is connected to the mains (for charging) or in operation.
- Never charge the device in a flammable environment (e.g. flammable objects, gasoline, flammable gases, paint, etc.)
- Only connect the unit for charging so that the cable and power connector are freely accessible.
- The built-in batteries are not intended for replacement or removal by the end user. If problems occur with the power supply, contact the customer service.
- This device must only be supplied at safety extra low voltage (SELV) corresponding to the marking on the appliance.
- To disconnect the unit from the mains, always pull the plug out of the power connection. Do not pull on the cable!

Before using the unit

- Check the device carefully for damage before each use. Do not operate a defective device.
- Do not use the device if any damage is visible, if it is not working properly, if it has been dropped or if it has become wet. To avoid hazards, send the unit to the service centre for repair.
- This product has strong light, please purchase qualified goggles, protective glasses can effectively protect 620nm wavelength of strong light, its light transmittance of 10%.

Intended Use

The **medisana** Radiance Trio LED Facial Device **RF 300** is a portable home beauty instrument designed for skincare. Enjoy spa-like treatments in the comfort of your home to enhance your natural beauty effortlessly.

- **RF, EMS and Gentle red-light:** Embrace skin radiance

Items supplied and packaging

Please check first of all that the device is complete and is not damaged in any way.

In case of doubt, do not use it and contact your supplier or your service centre.

The following parts are included:

- 1 **medisana** Radiance Trio LED Facial Device **RF 300**
- 1 USB-C charging cable
- 1 instruction manual
- 1 charging stand

Charging the battery

- Use a 5V main adapter and connect the supplied USB charging cable to the USB-C port of the charging stand **11**. Put the device on the charging stand and make sure the connector **8** on the device is aimed at the charging port **10** of the charging stand.
- During charging, the charging indicator **7** on the device will be lighted up in red. The charging time is about 2.5 hours. The device cannot be used during the charging process. When the battery is fully charged, the charging indicator **7** will turn blue.
- The device can also be charged directly without the stand by connecting the USB-C cable to the USB-C charging port **9** of the device.

Using the device

- Clean and dry your skin before using the device. To avoid stinging, wipe moisture off your face.
- Press and hold the Power/Mode button **6** for 2 seconds to switch on the device. The device starts working at RF mode. The RF mode indicator **5** will be lighted up in blue.
- Press the Power/Mode button **6** again to enter RF+EMS mode. The EMS mode indicator **5** will be lighted up in blue. During RF+EMS mode, press the EMS level button **4** to switch the EMS intensity level between Level 1, 2, 3. The corresponding EMS level indicator **3** will be lighted up in blue.
- After selecting the mode, apply the RF and EMS head **2** to the skin and slowly slide it. After 5 minutes of constant use, the device will stop and switch off automatically.
- Apply skin moisturizer on your face after the treatment.



CAUTION

Do not store or use the device in places with high temperature and humidity such as sauna and bathroom, etc.



WARNING

- **When the device comes into contact with the skin, red light 1 will be activated. Please do not use this product on the eyes!**
- **Do not stay in the same place for too long, RF itself will be hot. If stay too long, it is easy to cause skin burns**
- **If the device does not stop automatically after 5 minutes of use, do not use it and contact after-sales treatment.**

Troubleshooting

- Skin becomes reddened or trouble occurs after using the product.
 - Sensitive skin may experience some problems after using the product which may only be temporary. It is recommended to discontinue use for 2-3 days before use. If the redness and swelling do not improve after sufficient rest, please consult a specialist before using.
- Can the device also be used while charging?
 - It will be locked automatically in the state of charging and can not be used.
- I experience warmth while using the product.
 - It is completely normal to experience warmth while using, there may be difference in the temperature you feel depending on the skin type and sensitivity. However, if you experience inconvenience due to excessive heat, please stop using the product and seek consultation from the customer service center.

Maintenance and cleaning

- Children shall not play with the appliance and cleaning. User maintenance shall not be made by children without supervision.
- Do not store the product in places with high humidity, high temperature, a lot of dust or direct sunlight.
- After each use, use a dry, clean cloth and a specially formulated cleaner for electronic equipment to gently wipe the outer shell and the skin contact surface of the device.
- Do not store in places with presence of magnetic or electromagnetic wave.
- Do not drop the main body from high elevation or exert strong impact to the main body. If such occurs, make sure to have the product inspected before using the product.
- Do not disassemble the main body or attempt to repair at home. If there is breakdown, stop using the product and contact the store where you purchased the product.
- Dropping or exerting impact may become the cause of damage and breakdown since precision components have been built in the product.
- Since there is generation of high voltage within the main body, do not disassemble the main body arbitrarily.
- The LED of this appliance is not replaceable.

Disposal



This product must not be disposed of together with domestic waste.

All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your local authority or your supplier for information about disposal.

Technical specifications

Name and model:	medisana Radiance Trio LED Facial Device RF 300
Power supply:	Li-Ion battery 3.7V, 850mAh Input: 5V == 0.4A
Max. operating time:	5 minutes, with automatic switch-off
Operating conditions:	+10°C to +35°C, 3% to 75% max. relative humidity
Storage conditions:	-40°C to +70°C, 10% to 90% max. relative humidity
Dimensions (L×W×H):	Approx. 18.7 x 5.2 x 3.5 cm
Weight:	Approx. 115 g
Article No.:	50028
EAN No.:	4897138 50028 3

In accordance with our policy of continual product improvement, we reserve the right for technical and design changes without notice.

The current version of these instructions for use can be found at www.medisana.asia

Warranty and repair terms

Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is. The following warranty terms apply:

1. The warranty period for **medisana** products is one year from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

