


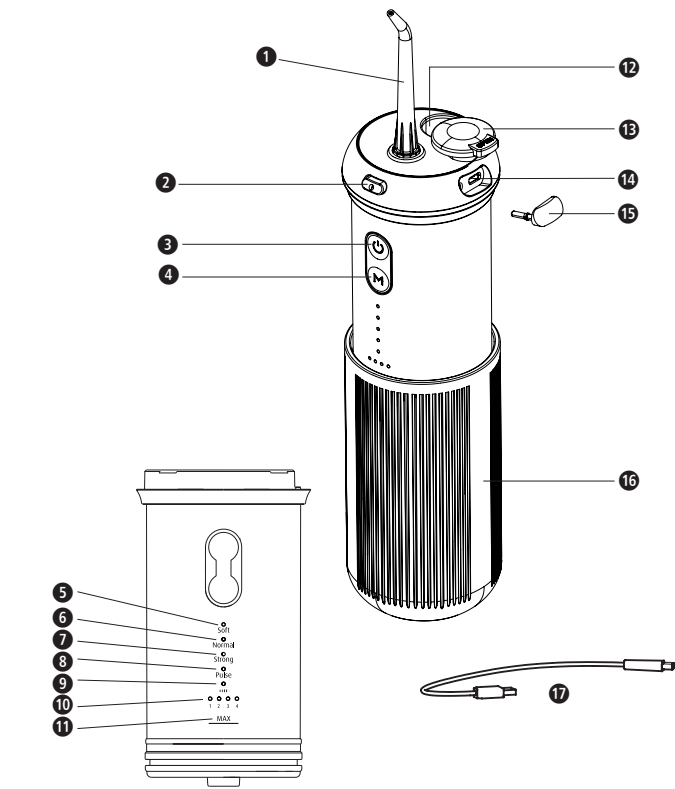

Device and controls

- 1 Nozzle
- 2 Nozzle release button
- 3 Power On/Off button
- 4 Mode button
- 5 Soft mode indicator
- 6 Normal mode indicator
- 7 Strong mode indicator
- 8 Pulse mode indicator
- 9 Battery indicator
- 10 Water pressure indicator
- 11 Maximum expansion point
- 12 Water injection port and nozzle storage hole
- 13 Top cover
- 14 USB-C charging port
- 15 Charging port cover
- 16 Water tank
- 17 USB-C charging cable

Explanation of symbols

This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.

- WARNING**
These warnings must be observed to prevent possible injury to the user.
- ATTENTION**
These instructions must be followed to prevent possible damage to the appliance.
- NOTE**
These notes provide you with useful additional information on installation or operation.

LOT LOT number

Recycling symbols/codes:
These are used to provide information about the material and its proper use and recycling.

Scope of delivery and packaging

Please check first of all that the device is complete and is not damaged in any way. In case of doubt, do not use the appliance and contact your dealer or your service centre. The following parts are included:

- 1 **medisana** Portable Water Flosser **WF 150**
- 1 Charging cable
- 2 Nozzles
- 1 Instruction manual

The packaging can be reused or recycled. Please dispose properly of any packaging material no longer required. If you notice any transport damage during unpacking, please contact your dealer without delay.

WARNING
Make sure that the packaging film does not get into the hands of children.
There is a risk of suffocation!

GB



IMPORTANT INFORMATION! RETAIN FOR FUTURE USE!

Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. Should you give this device to another person, it is vital that you also pass on these instructions for use.



Safety Information

- When using an electronic appliance (especially if children are present), basic precautions should always be followed, including the following.
- Use the device only according to its intended purpose as specified in the instruction manual. The warranty will be invalidated if used for purposes other than those for which it is intended.
- Any application other than the herewith described is improper and therefore considered dangerous. If the appliance is not used according to the instructions, the user is liable for the safe operation of the unit.
- Water flosser can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- The use of the appliance on children under 8 years of age is not recommended. Children aged 8 and over must be supervised by an adult.
- It is also recommended that adults first learn the standardised use of water flosser before instructing children in their use.
- We recommend that children only ever use the soft intensity mode and the lowest pressure level.

Caution! Reduction of the risk of electric shock:

1. Do not use when bathing.
2. Never store the appliance in such a way that it can fall or be pulled into a bath, shower cubicle or sink.
3. Do not place or drop in water or other liquids.
4. Do not reach for an appliance that has fallen into water. Unplug the appliance immediately.
5. Do not touch the charger with wet hands when charging.

Warning! To reduce the risk of burns, electric shock, fire or injury:

1. The appliance should never be left unattended.
2. Close supervision is required when this product is used by, on or near children or disabled persons.
3. Use the appliance only for the purpose described in this manual. Do not use attachments that are not recommended by the manufacturer.
4. Do not operate the appliance if it is damaged, if it is not working properly, if it has been dropped or damaged or if it has fallen into water. Take the appliance to a service centre for inspection and repair.
5. Keep the cable away from hot surfaces.
6. Never block the air openings of the appliance and do not place it on a soft surface, such as a bed or couch, where the air openings could be blocked. Keep the air openings free of fluff, hair and the like.
7. Never use the appliance when you are asleep or drowsy.
8. Do not drop or insert any objects into the openings.
9. Disconnect the charging cable before filling the appliance. Only fill the device (water tank 16) with water, unless otherwise specified by the manufacturer. Do not overfill!
10. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
11. Children should be supervised to ensure that they do not play with the appliance.
12. Do not use the device outdoors or in places where aerosol products (sprays) are used or where oxygen is administered.
13. The appliance may only be operated with safety extra-low voltage in accordance with the labelling on the appliance.
14. Do not connect this device to a voltage system other than that specified on the device or charger.
15. This device contains batteries (rechargeable battery) which are not replaceable.
16. The device is only suitable for domestic use. Do NOT use the product outdoors.

Warning! To reduce the risk of electric shock and other hazards, please observe the following safety guidelines:

- Check whether the charger or the cable is damaged before you start up or continue to use the device. If either is damaged, do not attempt to use the device. Instead, contact customer service for a replacement.
- Do not throw the appliance into a fire or other hot environment.
- Do not immerse the device in water as this may cause damage.
- Do not attempt to plug or unplug the charger with wet hands.
- Do not modify, damage or destroy the device. Do not pull, bend or twist the cable by force. Do not place heavy objects on the cable and do not clamp the cable between other objects.
- Do not disassemble the device, except to remove the battery for disposal (see corresponding chapter).
- Do not place the device, the charger or the cable on heated surfaces.
- Do not place the device in a hot or humid place or in a place exposed to direct sunlight.
- Before cleaning and maintenance, ensure that the device is switched off and is not charging.

Usage: Charging the battery


We recommend fully charging the **medisana WF 150** Water Flosser before using it for the first time (maximum charging duration: approx. 3 hours). The device cannot be used during the charging process. To charge, connect the supplied cable to the charging port 14 and to a 5V /1A power source. During the charging process, the battery indicator 9 flashes. When the battery is fully charged, the battery indicator 9 lights up green..

If the battery level is low during operation, the battery indicator 9 flashes slowly - charge the battery as soon as possible. If the battery charge level is so low that further use is not possible, the battery indicator 9 flashes in quick succession before the appliance is automatically switched off. The appliance cannot be used again until the battery has been recharged.

Installing / Removing the spray nozzle

Insert the nozzle into the centre hole of the device unit it clicks into place.

When the nozzle needs to be removed or replaced, it can be easily removed just by pressing the nozzle release button 2 and pulling out the nozzle at the same time.

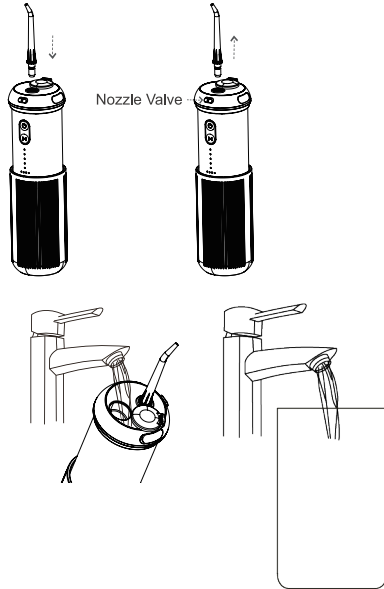


Filling the water tank 16

Open the top cover 13 and expand the device by pulling the water tank 16 down from the device until reaching the maximum expansion point 11. Fill water into the water injection port 12. Alternatively, you can carefully detach the water tank from the device (see drawing below) to fill it more quickly.

Note: Make sure that the device is switched off when filling water tank 16.

Then reconnect the water tank to the device.

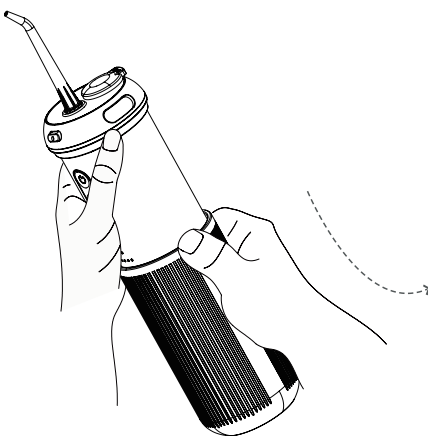

Detach the water tank 16 from the device:
Open the top cover 13 before pulling the water tank 16 so that it slides smoothly. Then, grasp the device housing with your left hand and press the water tank slightly downwards with your right index finger; at the same time, pull firmly anti-clockwise downwards to the right to remove the water tank.

Operation

Power On/Off button 3 :

- Press the Power On/Off button 3 to switch on the device. Press the button again to switch the device off.

Mode button 4 :


- Press the mode button 4 to switch between 4 intensity modes: Soft, Normal, Strong, Pulse.

- Each intensity mode includes 4 adjustable pressure level 10. Press and hold the mode button 4 for 2 seconds to switch the water pressure level 10 of current mode.

NOTE: If you are using the device for the first time, it is recommended to start at the soft mode. You can increase the level later depending on your personal preference.


Using Water Flosser

1. Lean over the wash basin and put the nozzle tip 1 in the mouth. Aim the tip at your teeth.

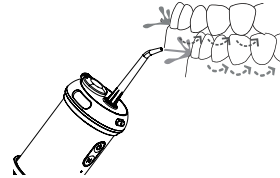

2. Turn on the device by pressing the power On/Off button 3. Press the mode button 4 to switch to the suitable mode and pressure level.

3. Rinse the gingiva at 90° and close your lips slightly to avoid water from splashing out of the mouth. For best results, start with the molars and clean towards the front teeth. Slide the nozzle tip along the gumline and keep flossing a while at interdental space.

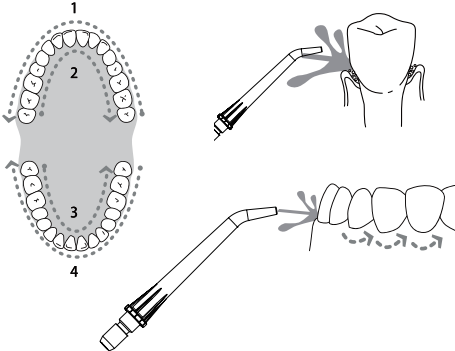

After use

Switch off the appliance, pull the water tank 16 down until it is separated from the main body. Pour out the remaining liquid and install the water tank back to the device.

Cleaning and care

Cleaning the appliance:

- Only clean the device with clean water or a neutral cleaning agent. Do not use any corrosive or abrasive cleaning agents.

- Never place the device in water or liquids. If necessary, clean the surface of the device with a soft cloth dipped in a mild, non-corrosive cleaning agent.

Cleaning the water tank 16:

- Hold the device and the water tank 16 with your right hand. Pull the water tank 16 downwards to separate it from the device.
- Rinse the inside of the water tank 16 with clean water at room temperature. Do not clean the water tank 16 with additional cleaning agents.
- After cleaning, put the water tank 16 back in its original position.

Removal of mineral deposits caused by hard water:

Depending on the mineral content of the water used, deposits may form on the device. If you do not remove deposits from hard water, the performance of the device will be impaired. We therefore recommend repeating the following cleaning process every 1-3 months to maintain the performance of the appliance:

Point the spray nozzle at a sink, switch on the device and let it run until all the water has been used up. Fill water tank 16 with clean, hot water (>85°C) and then repeat the above procedure.

Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Please remove the batteries before disposing of the device/unit. Do not dispose of old batteries with your household waste, but at a battery collection station at a recycling site or in a shop. Consult your municipal authority or your dealer for information about disposal.

Troubleshooting

Problem	Possible cause	Solution
Device is leaking (e.g. leak between spray nozzle and device body or in water tank 16)	The seal is damaged or the spray nozzle has not been installed in the correct position.	Install the spray nozzle or the water tank 16 in the correct position.
Insufficient spray pressure	Battery not sufficiently charged.	Charge the battery.
Device does not work	Battery not sufficiently charged.	Charge the battery.
Device is working, but no liquid comes out of the spray nozzle	Device is being used for the first time or has been used for a long time.	Fill the water tank 16 to the maximum and allow the appliance to operate at the highest level for a few seconds.

If the error cannot be rectified using the table, please contact customer service.

Technical specifications

Name and model:	medisana Portable Water Flosser WF 150
Power supply:	3.7V== 1000mAh / Input: 5V== / 1A 5W (via USB)
Power consumption:	10W
Noise level:	< 70 dB
Water tank capacity:	190 ml
Water resistance:	IPX7
Charging time:	Approx. 3 hours
Operating time:	30 minutes for Soft mode (Auto shutdown in 2 minutes, can be used 15 times in 1 charge)
Dimensions approx.:	Approx. 58 x 58 x 133 mm
Weight (main unit):	Approx. 230 g
Article no.:	50159
EAN code:	48 97138 50159 4

In accordance with our policy of continual product improvement, we reserve the right to make technical and optical changes without notice.

Warranty and repair terms

Please contact your supplier or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is. The following warranty terms apply:

1. The warranty period for **medisana** products is 1 year from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
2. Defects in material or workmanship will be removed free of charge within the warranty period.
3. Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
4. The following is excluded under the warranty:
 - a. All damage which has arisen due to improper treatment, e.g. nonobservance of the user instructions.
 - b. All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - c. Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - d. Accessories which are subject to normal wear and tear.
5. Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

medisana.

MEDISANA MALAYSIA SDN BHD

G-02, H2O Residences, Jalan Pju 1a/3, Ara Damansara, 47301 Petaling Jaya, Selangor Darul Ehsan, Malaysia